# T4 VET control kit

**For veterinary use only!**

For quality control procedures of the quantitative determination with the T4 VET test kit on a solo or CUBE-VET analyser

Eurolyser Diagnostica GmbH  
Bindergasse 3  
5020 Salzburg/Austria

Tel. +43 (0)662 / 43 21 00  
Fax +43 (0)662 / 43 21 00-50  
www.eurolyser.com

## English

|   | Order information   | Indication   | Kit size   |
|---|---|--|--|
|  | Order number: VT 0140<br>Order number: VT 0141<br>Order number: VT 1400 | T4 VET test kit<br>T4 VET test kit<br>T4 VET control kit | 16 tests/pack<br>6 tests/pack<br>1 x 5 ml (decision level) |

**Test kit preparation: Allow single test at least 10 minutes to warm up to room temperature (20 - 25 °C) by placing the test into the test kit rack. Put test kit package back into refrigerator.**

### Summary

This T4 control kit is intended for use as a test sample of known concentration for the solo T4 assay. The user can compare observations with expected ranges as a means of assuring consistent performance of the solo T4 assay. The control is prepared from a bovine serum base. Analyte levels are adjusted with various drugs, drug metabolites and purified chemicals. Preservatives and stabilisers are added to maintain product integrity.

### Required Material

*(not supplied with the kit)*

- Pipettes

### Preparation and Handling

The control is liquid and ready to use. Thoroughly mix the contents of the vial before each use by gently inverting for several minutes. Open the vial and transfer the required quantity of control into a clean sample cup. Use sample cup only for one measurement circle and use new sample cups further on. Replace the cap immediately and store the opened vial at 2 - 8 °C. Always use disposable pipette tips!

### Stability and Storage

Unopened vials are stable at 2 - 8 °C until the expiration date stated on the label. Once opened, control vials are stable for 30 days when stored tightly capped at 2 - 8 °C. Bacterial contamination produces an increase in turbidity and/or a characteristic odor. Discard vial if evidence of microbial contamination is observed. **DO NOT FREEZE!**

### Waste Management

Please refer to local legal requirements

### Warnings and Precautions

This control kit is for veterinary diagnostic use only!  
**DO NOT INGEST!** Exercise the normal precautions required for handling all laboratory reagents. Components of the control which are derived from human source material have been tested using FDA accepted methods and found nonreactive for Hepatitis B Surface Antigen (HBsAg), Hepatitis C (HCV), HIV1 and HIV2. However, no test method can offer complete assurance that products derived from human source material are free of infectious agents.

Set the species to „Control“ before running the QC material

Compare the results with the values indicated in the lot-specific value table

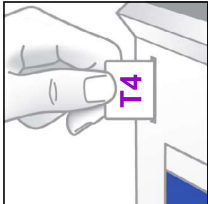
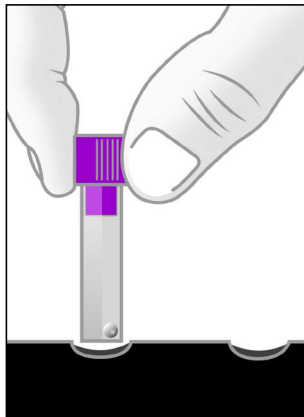
If the result is not within the acceptable range, repeat measurement. If result is again not within the acceptable range, please contact your distributor

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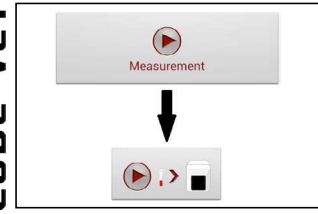
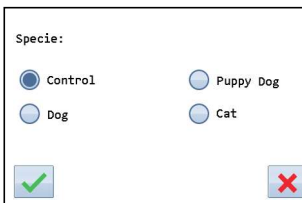
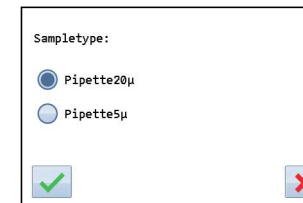
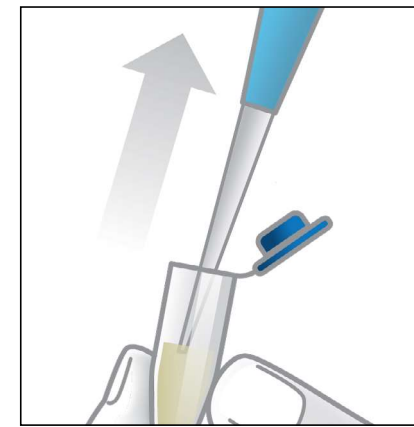
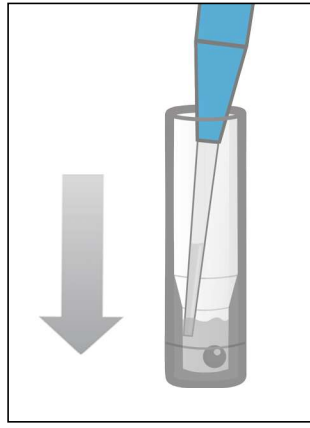
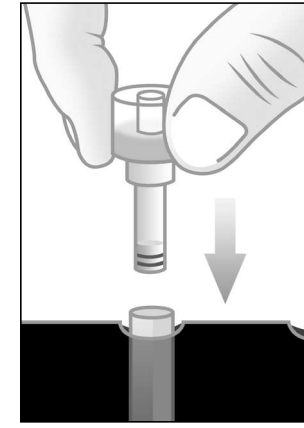
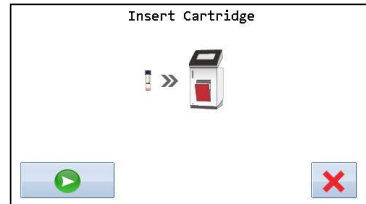
## English

### ATTENTION!

Allow single test at least 10 minutes to warm up to room temperature (20 - 25 °C) before use!

- 1. Preparation of test system**
  - 1.1 Place RFID card
  - 1.2 Place ERS cuvette in test kit rack
  - 1.3 Place ERS cap in test kit rack
  - 1.4 Press „Measurement“ button, enter required information using the touchscreen

- 2. Control preparation**

*For further details please see user manual of laboratory photometer*

  - 2.1 Set species to „Control“
  - 2.2 Set sample type to „Pipette 20µ“
  - 2.3 Aspirate 20 µl control liquid using a pipette

- 3. Control processing**
  - 3.1 Dispense 20 µl control liquid INTO THE LIQUID in the ERS cuvette
  - 3.2 Apply ERS cap firmly onto ERS cuvette
  - 3.3 Place ERS cartridge into laboratory photometer
  - 3.4 Start automatic control processing by pressing the  start button on the solo laboratory photometer, or by closing the door of the CUBE-VET laboratory photometer.

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Technical details subject to change without notice.

Doc-ID: d68795b0-2f0e-448a-8008-f810ddb79672 ; latest revision: DP 2.0 ; date of print: 2018-06-13

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